







Priorclave's Base range covers "the basics": media preparation, waste loads, and sterilising loose instruments, open containers, and labware. These are high-efficiency cylindrical chambered autoclaves, caster-mounted for easy transport. They are well suited to research, quality assurance, and teaching labs. Easy maintenance and light resource consumption help ensure a low lifetime cost of ownership. All Base models include fully programmable cycle customisation, "media warming," integrated cooling and freesteaming systems, and drain protection.

SPEAK WITH AN AUTOCLAVE EXPERT







+44 (0)20 8316 6620 • priorclave.co.uk • sales@priorclave.co.uk


Our highly trained service technicians provide dedicated pre- and after-sales support including planned maintenance programmes to ensure you get the most from your Priorclave.

KEY FEATURES

- Quickseal single-action door fitted with thermal and pressure locks to prevent opening at load temperature above 80°C and pressures above 0.2 bar.
- Tactrol® 3 microprocessor control with simple and fully variable setting of process time, temperature with graphic indication of cycle status.
- · Automatic timed free-steaming.
- Media Warming and Delayed Start functions keep autoclaved media at a usable temperature ready for immediate use.
- · Low water level sensor.
- Two full width stainless steel load shelves.
- Epoxy coated panels and frame members are treated with an anti-bacterial agent which is effective against all bacteria and fungi including MRSA.
- Pressure vessel carries a 10 year warranty and is insurance approved.
- Service and support direct from manufacturer (or appointed distributor).

320L AUTOCLAVE FRONT LOADING


EH320-BASE, SH320-BASE

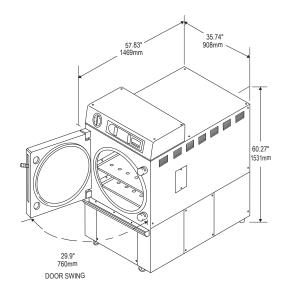

Heat Input	Chamber Volume	Operating Range	Loading Format	Loading Height	Heating Power	
Internal Electrical or Direct Steam Source	320 litres	Up to 138°C, 2.4 Bar (280°F 34 psi)	Front Loading	725mm (X")	18kW 30A three-phase	
Chamber Dimensions diameter x depth			Bottom Shelf Dimensions w x d		Top Shelf Dimensions w x d	
630 x 1025mm (24.80" x 40.35")		130 %	490 x 1025mm (19.29"x 40.35")		620 x 1025mm (24.40" x 40.35")	
		Weight (unloaded)	Minimum Install	ation [Door Swing	
908 x 1469 x 1531mm (35.74" x 57.83" x 60.27") (300 kg (661.38 lbs)	Please allow 500mm (either side and 300mm at the back of the Aut	(11.81")	760mm (29.9")	

400L AUTOCLAVE FRONT LOADING

EH400-BASE, SH400-BASE

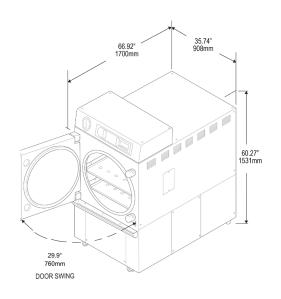

Heat Input	Chamber Volume	Operating Range	Loading Format		ding ght	Heating Power
Electrical or Direct Steam Source	400 litres	Up to 138°C, 2.4 Bar (280°F 34 psi)	Front Loading	725	mm	18kW at 30A three phase
Chamber Dimensions diameter x depth		Bottom Shelf Dimensions w x d		Top Shelf Dimensions w x d		
630 x 1200mm (13.8" x 24.4")		490 x 1200mm (19.29" x 47.24")		620 x 1200mm (24.40" x 47.24")		
Outside Dimensions w x d x h		Weight (unloaded)	Minimum Installation		Door Swing	
908 x 1700 x 15: (35.74" x 66.92" x	J	350kg (771.61 lbs)				760mm (29.9")

CONFORMITY & COMPLIANCE

UK Standard	ISO 9001 / CE Marked / BS2646
US Standard	ASME
EU Standard	CE / ISO 9001 / PED 97/23/EC

INT Standard	EMC Compliance
Sustainability	ACT Label Issued
Service	Priorclave Trained Standard