



Priorclave's Base range covers "the basics": media preparation, waste loads, and sterilising loose instruments, open containers, and labware. These are mid-sized, efficient rectangular chambered autoclaves designed for research, quality assurance, and teaching labs. All Base models include fully programmable cycle customisation, "media warming," integrated cooling and freesteaming systems, and drain protection.

KEY FEATURES

- Lower handwheel on door is fitted with thermal and pressure locks to prevent opening at load temperature above 80°C and pressures above 0.2 bar.
- Tactrol® 3 microprocessor control with simple and fully variable setting of process time, temperature with graphic indication of cycle status.
- Automatic timed free-steaming.
- Low water level sensor.
- Two full width non-tipping stainless steel load shelves.
- Epoxy coated panels and frame members are treated with an anti-bacterial agent which is effective against all bacteria and fungi including MRSA.
- Pressure vessel carries a 10 year warranty and is insurance approved.
- Service and support direct from manufacturer (or appointed distributor).

SPEAK WITH AN AUTOCLAVE EXPERT




+44 (0)20 8316 6620 • priorclave.co.uk • sales@priorclave.co.uk

Our highly trained service technicians provide dedicated pre- and after-sales support including planned maintenance programmes to ensure you get the most from your Priorclave.

230L AUTOCLAVE RECTANGULAR


EH230-BASE, SH230-BASE

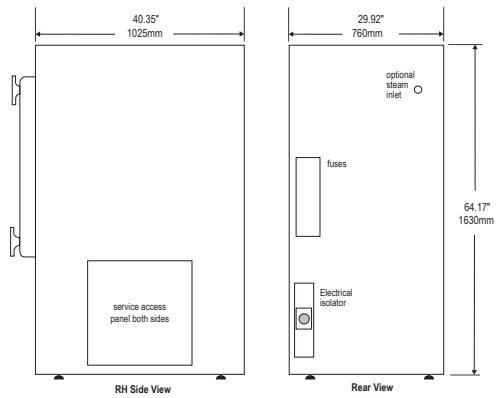

Heat Input	Chamber Volume	Operating Range	Loading Format	Loading Height	Heating Power
Electrical or Direct Steam Source	230 litres	Up to 138 °C, 2.4 Bar (280 °F 34 psi)	Front Loading	800mm (31.49")	21kW 30A three-phase
Chamber Dimensions w x d x h		Bottom Shelf Dimensions w x d		Top Shelf Dimensions w x d	
580 x 600 x 600mm (22.83" x 23.62" x 23.62")		565 x 550mm (22.24" x 21.65")		565 x 550mm (22.24" x 21.95")	
Outside Dimensions w x d x h	Weight (unloaded)	Minimum Installation w x d x h	Door Swing		
760 x 1025 x 1630mm (29.92" x 40.35" x 64.17")	560 kg (1234.59 lbs)	1225 x 960 x 760mm (48.22" x 37.79" x 29.92")	760mm (24.6")		

285L AUTOCLAVE RECTANGULAR


EH285-BASE, SH285-BASE

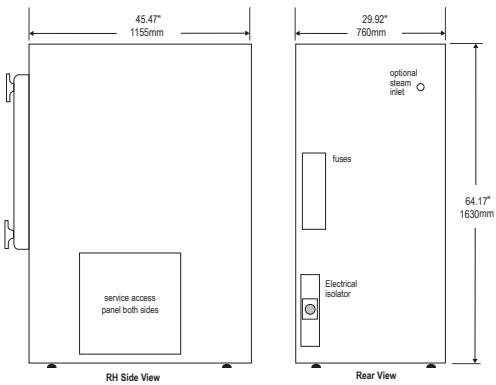

Heat Input	Chamber Volume	Operating Range	Loading Format	Loading Height	Heating Power
Electrical or Direct Steam Source	285 litres	Up to 138 °C, 2.4 Bar (280 °F 34 psi)	Front Loading	800mm (33.9")	21kW 30A three-phase
Chamber Dimensions w x d x h		Bottom Shelf Dimensions w x d		Top Shelf Dimensions w x d	
580 x 760 x 600mm (22.83" x 29.92" x 23.62")		565 x 660mm (22.64" x 25.98")		565 x 660mm (22.24" x 25.98")	
Outside Dimensions w x d x h	Weight (unloaded)	Minimum Installation w x d x h	Door Swing		
760 x 1155 x 1630mm (29.92" x 45.47" x 64.17")	600 kg (1322.77 lbs)	1355 x 960 x 760mm (53.34" x 37.79" x 29.92")	760mm (24.6")		

350L AUTOCLAVE RECTANGULAR

EH350-BASE, SH350-BASE

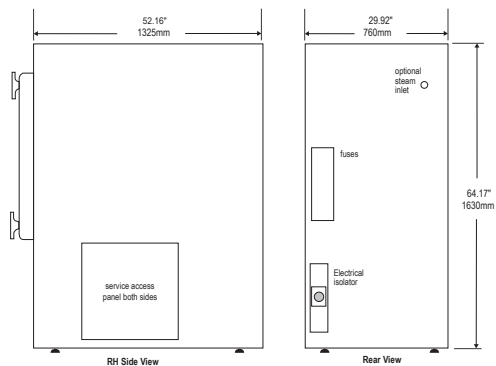

Heat Input	Chamber Volume	Operating Range	Loading Format	Loading Height	Heating Power
Electrical or Direct Steam Source	350 litres	Up to 138 °C, 2.4 Bar (280 °F 34 psi)	Front Loading	800mm (33.9")	21kW 30A three-phase
Chamber Dimensions w x d x h		Bottom Shelf Dimensions w x d		Top Shelf Dimensions w x d	
580 x 900 x 600mm (22.93" x 35.43" x 23.62")		565 x 850mm (22.24" x 33.46")		565 x 850mm (22.24" x 33.46")	
Outside Dimensions w x d x h	Weight (unloaded)	Minimum Installation w x d x h	Door Swing		
760 x 1325 x 1630mm (29.92" x 52.16" x 64.17")	640 kg (1410.96 lbs)	1525 x 960 x 760mm (60.03" x 37.79" x 29.92")	760mm (24.6")		

CONFORMITY & COMPLIANCE

UK Standard	ISO 9001 / CE Marked / BS2646
US Standard	
EU Standard	CE / ISO 9001 / PED 97/23/EC

INT Standard	EMC Compliance
Sustainability	
Service	Priorclave Trained Standard