


Autoclave Design Differences

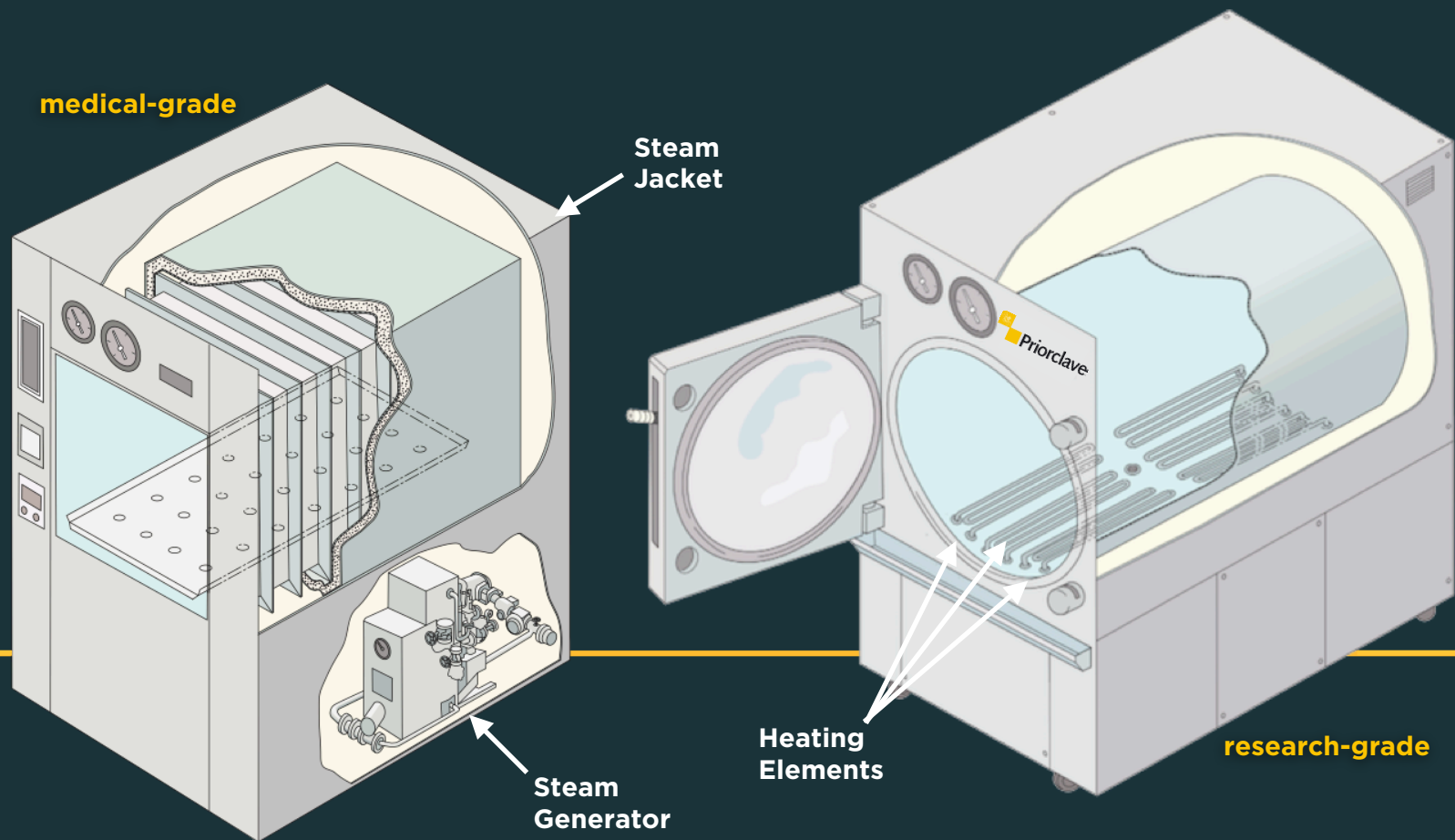

Rectangular Model with Steam Jacket

- Steam is generated by external boiler
- Reinforced corners in chamber
- Great for high throughput (fast!)
- Designed for sterilization (flat instruments)
- Requires RO/DI/Hot water connection

Cylindrical Model (No Steam Jacket required)

- Steam is generated by elements inside chamber
- Natural steam circulation (no cold pockets)
- Ideal for up to 5 cycles/day (on-demand)
- Designed for media, glassware, waste
- Requires soft tap water connection or manual soft water feed