

## Autoclave Consumables

Available direct from Priorclave, one of the UK's top autoclave manufacturers, is a comprehensive selection of competitively priced consumables all designed to help you maintain your steam steriliser in tip-top condition and ensure perfect operation.

### Syphon Pumps

This range of hand syphon pumps available in two sizes is designed to make light work of removing water from the bottom of top and front loading autoclaves.

Description	Part number
Short Syphon Pump for front loading autoclaves	ZZZ/SYP/SML
Long Syphon Pump for top loading autoclaves	ZZZ/SYP/LGE



### Deodorising Capsules

Priorclave deodorising capsules containing concentrated perfumed essential oils come in a choice of three different fragrances (packs of 100).

Description	Part number
Classic - Deodorising Capsules	ZZZ/DEO/001
Apple - Deodorising Capsules	ZZZ/DEO/002
Citrus - Deodorising Capsules	ZZZ/DEO/003



### Printer Ribbons and Rolls

Ensure smart legible printouts with Priorclave printer ribbon replacement units.

Description	Part number
Type 1 Printer Ribbon (Straight)	ZZZ/RIB/001
Type 2 Printer Ribbon (L Shaped)	ZZZ/RIB/002
Printer Paper Rolls	ZZZ/ROL/005
Thermal Printer Rolls	ZZZ/ROL/007

type 1



type 2



### Bowie & Dick Test Pack

It is recommended that a steam penetration test is run in every autoclave; with the Bowie & Dick test the colour of the indicator sheet changes from Blue to Pink confirming that steam penetration has been effective up to the centre of the sterilisation pack.

Description	Part number
BD-1XE Bowie & Dick Test (20 sheets per pack)	ZZZ/BDT/000



Priorclave Limited  
 129-131 Nathan Way, West Thamesmead Business Park  
 London, SE28 0AB  
 T: +44 (0)20 8316 6620 F: +44 (0)20 8855 0616  
 E: service@priorclave.co.uk

[www.priorclave.co.uk](http://www.priorclave.co.uk)

## Biological Indicators

Now available from Priorclave are Self-Contained Biological Indicators (SCBI) and Mini Self-Contained Biological Indicators (MSCBI) for monitoring steam sterilisation processes (spore population = 10<sup>5</sup> and 10<sup>6</sup>).

Each SCBI unit consists of a plastic container with a cap, a crushable glass ampoule with recovery media and a disc inoculated with spores. They are certified for species, population, purity, resistance (d value, z value, where applicable) and expiration date.

Description	Part number
Self contained Geobacillus stearothermophilus (boxes of 50)	ZZZ/SCS/06E
Mini self contained Geobacillus stearothermophilus (boxes of 100)	ZZZ/MSS/06E



## Spore Ampoules

Spore ampoules are used for monitoring the efficacy and validation of the steam sterilization processes. No activation is required; simply expose spore ampoule to sterilization process and incubate. All spore ampoules, available in several configurations, are hermetically sealed type I borosilicate glass containing a modified Soy-bean Casein Digest Broth with bromocresol purple pH indicator and the specified population of Geobacillus stearothermophilus spores. Growth is evident by either turbidity and/or a colour shift from a purple to or toward a yellow colour.


Description	Part number
Spore Ampoules (box of 15)	ZZZ/SA1/05E



## TST Strips

The TST (Time, Steam and Temperature) steriliser control is the most sophisticated chemical device that indicates a PASS only when it has been subject to the required heat and time ratio in saturated heat steam when exposed to steam for 3-5 minutes at 134°C or 15 minutes at 121°C.

Description	Part number
TST Strips (boxes of 100)	ZZZ/TST/S-E

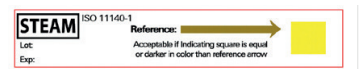

## Adhesive tape

White adhesive tape for use with standard autoclaves. It reveals the word 'Dark is Exposed' after 15 minutes at 121°C and provides clean removal from the item. Tape sizes: 15mm x 50m, blue to dark grey/black.

Description	Part number
Adhesive tape	ZZZ/TAP/S-E



Priorclave Limited  
129-131 Nathan Way, West Thamesmead Business Park  
London, SE28 0AB  
T: +44 (0)20 8316 6620 F: +44 (0)20 8855 0616  
E: service@priorclave.co.uk

[www.priorclave.co.uk](http://www.priorclave.co.uk)